

# Electronic pipette controller

## The Profiller™ *electro* 447

**This pipette controller fits any hand in a very comfortable way. Shape is ideal for working with glass or plastic pipettes within laminar flow or microbiological safety cabinet. A maximal working efficiency, even with large pipetting volumes, is reached through a powerful, yet quiet motor pump and the availability of several setting options. One-year warranty.**

- Powerful pump – fills 25 mL in 3 seconds
- Comfortable finger triggers
- Large display always visible
- Real time battery charge indicator
- NiMH battery rechargeable during use
- Reliable construction for shock and chemical resistance
- Interchangeable protection membrane filter
- Compact table and wall stands



**447**

Pipettes  
up to 100 mL and over



### Choice of dispensing modes

Instant activation of dispensing with blow-out (TC) or gravity mode (TD).


### Instant speed setting

A switch enables selection of fast or slow aspiration and dispensing speeds, clearly indicated on display.

### Ordering information

Profiller™ <i>electro</i> 447	Packaging	Cat. No.
With plug and charger	1 / pk	447:100

Description	Packaging	Cat. No.
Membrane filter, 0.45 µm	5 / pk	322.447
Membrane filter, 0.2 µm	5 / pk	322.407
Silicone pipette holder	1 / pk	1.447.01
Nose cone	1 / pk	1.447.02

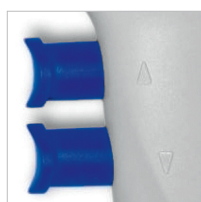

### Finger triggers

Aspiration and dispensing controlled through two ergonomically positioned finger triggers requiring minimum activation effort. Triggers enable gradual speed regulation.



### Battery charge

Instrument plug enables charging battery either in storage position or during use. Charge level visible on display.

### Convenient stands

Instrument supplied with stable bench stand including 3 suckers, and adhesive wall hanger.

### Supply Content

Supplied with power supply and cable, four plugs adaptors (EU, USA/JP, GB, AUS/NZ), two spare hydrophobic filters (0.45 µm and 0.2 µm), bench stand, wall hanger and instruction manual.

Description	Packaging	Cat. No.
Spare bench stand	1 / pk	320.947
Spare wall hanger	1 / pk	320.447
Spare 1.2V NiMH battery	3 / pk	900.917