

Manual pipette controller


The Profiller™ manual 437

Robust and easy to use, this model offers ergonomic shape and lightweight. It provides smooth control in both aspirating and dispensing, using glass or plastic pipettes of volume up to 100 mL. One-year warranty.

- Large capacity bulb
- Soft action lever button
- Efficient blow-out system
- PTFE protection membrane filter
- Maintenance free

437

Pipettes
up to 100 mL

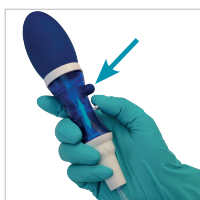

Squeezable bulb

Simply squeeze large capacity silicone bulb for efficient aspiration. If additional aspiration is needed, it can be squeezed even during pipetting.

Precise control

Thumb lever button for left and right-handed users to easily control both aspiration and dispensing speeds.

Liquid blow-out

Press large button to expel residual liquid when using blow-out type pipettes.

Protection filter

Interchangeable hydrophobic membrane filter guarantees effective protection against airborne contamination during aspiration and dispensing.

Pipette holder

Conical shape of silicone pipette holder tightly fits graduated glass or plastic pipettes. Autoclavable.

Ordering information - Profiller™ manual 437

Description	Packaging	Cat. No.
Profiller™ manual 437, blue bulb	1 / pk	437100

The Profillette™ 406

Low cost, easy to use pipette filler suitable for graduated glass and plastic pipettes. Rotate wheel to aspirate. Fast release by pressing air inlet valve button, or gradual release by gently turning rotation wheel.



406

2 mL
10 mL
25 mL

- One-hand operation
- Easy filling and dispensing
- Air inlet valve for fast distribution
- Colour coding per size
- Maintenance free

Ordering information - Profillette™ 406

Volume	Colour	Packaging	Cat. No.
2 mL	blue	1 / pk	406.002
10 mL	green	1 / pk	406.010
25 mL	red	1 / pk	406.025
2 / 10 / 25 mL	assorted	3 / pk	406.300

Ordering information - Accessories

Description	Packaging	Cat. No.
Profillette™ 406		
Spare silicone pipette holder	1 / pk	1.406.01
Replacement nose screw	1 / pk	1.406.02
Profiller™ manual 437		
Membrane filter, 0.45 µm	5 / pk	322.447
Silicone pipette holder	1 / pk	1.437.01
Nose cone - white	1 / pk	1.437.02